ThermoKlinik

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 Dear customer,
We thank you for having purchased a **Pic ThermoKlinik**, the thermometer that takes a fast, precise body temperature reading from a distance.

Pic ThermoKlinik is an innovative method for measuring body temperature. Thanks to the exclusive infrared technology it detects body temperature from a distance of 10 cm from the forehead.

Innovative technology, ultramodern design. ThermoKlinik is the Pic Indolor solution, ideal for the entire family.

All Pic Indolor thermometers are produced using the most advanced technologies and undergo rigorous quality controls, in compliance with European Community regulations, thus guaranteeing utmost quality and precision.

How it works

Using advanced infrared technology, Pic ThermoKlinik measures the temperature of the blood running through the temporal artery and converts this temperature into the equivalent of an oral temperature reading.

Characteristics

1. Probe
2. Display
3. Measurement Button
4. Memory Button
5. Protective, rotating cover

Body temperature

Like blood pressure, body temperature varies from person to person and undergoes different fluctuations during the course of the day, running approximately from 35.5°C to 37.8°C (95.5°F - 104.04°F).


For these reasons, to correctly determine the temperature, it is best to know one's normal, healthy forehead temperature at various times of the day. This will provide a valid reference when evaluating any real rise in temperature. To determine if one has a fever, compare the temperature detected with the person's normal temperature. Remember, a rise over the reference body temperature of 1° C (1.8°F) or more is generally indication of fever.

How to take a measurement

- 1) The thermometer must be set in an environment where the reading is to be taken 20-25 minutes before taking the reading. The person who will be taking the reading must also be in the room for at least 20-25 minutes.
- 2) Clean all sweat, make-up, etc. from the forehead and wait a few minutes before taking the reading.
- 3) Open the thermometer by turning the cover (fig. 1 - no. 5). The thermometer goes on automatically. A blue ray appears which helps ensure correct positioning of the probe when taking the measurement. In fact, this ray indicates the exact point where the reading will be taken. **Wait until the “READY” message appears before taking the reading.**

- 4) Position the thermometer's detector at a distance of not more than 10 cm from the forehead. Press the “Measurement” button (fig. 1 - no. 3) and explore the entire temple area. The best position is above the left side of the left eyebrow or above the right side of the right eyebrow, just under the hairline.

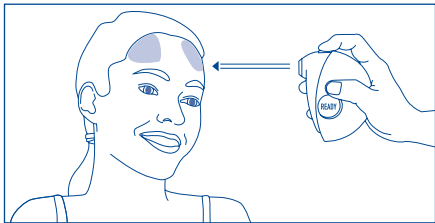While scanning, the thermometer emits short beeps to indicate that the temperature measurement is in progress.

- 5) The measurement is complete only after a long beep has sounded. At this point you can read the temperature on the display.
- 6) The thermometer turns off automatically 1 minute after the measurement has been completed. To turn off the blue positioning ray, turn the protective cover (remember, if the cover is not closed, the ray will continue to operate and use up the power supply batteries).

Note: Correct measurement can be taken at any distance up to 10 cm.

Important: For repeated measurements, always position the thermometer at the same point on the forehead. Wait at least 2 minutes between consecutive measurements and turn the thermometer off and on again.

Important notes

- After physical exercise, bathing or eating, wait 30 minutes before taking a temperature reading.
- Wait at least 2 minutes between consecutive measurements and the thermometer must be turned off and on again.
- When pointed toward a hot object, the thermometer detects the temperature remotely. To prevent incorrect readings, position the thermometer as soon as the “READY” message appears.
- Avoid direct sunlight, other sources of heat and direct air conditioning as they can affect the forehead during the temperature reading, modifying the result.
- If the light ray is no longer visible, replace its batteries.
- Never expose the thermometer to extreme temperatures, dust or abrupt movements.
- Never let the thermometer get wet.
- Never use the thermometer on irritated skin or when the forehead presents other problems.
- Never use the thermometer when the forehead is particularly sweaty, the reading will not be reliable.

Memory function

The memory function records 8 readings and remembers them for 64 hours. The last

reading performed before the thermometer is turned off is recorded to memory. To use the memory function, turn the protective cover (fig. 1 - no. 5), press the “M” button (located under the probe, fig. 1 - no. 4) and release it. The display alternates the following information four times:

- previous measurement
- how many hours ago the measurement was taken.

Pressing and releasing the “M” button again will call up the second to last measurement taken. You can read all the values saved by continuing to press the memory button.

To delete the memory, press the buttons “M” and “Measurement” simultaneously. This will turn off the thermometer and reset the memories.

Replacing the battery

Pic ThermoKlinik is powered by a 3V lithium battery (it also has two 1.5 V AAA batteries to power the light ray). The lithium battery can handle approximately 12,000 measurements. If the  symbol appears on the display or if the thermometer does not completely function, immediately remove the battery as indicated below:

- remove the lid from the battery compartment;
- remove the old battery;
- insert the new battery with the positive side (+) facing up and the negative side (-) facing down;
- replace the battery compartment lid.

Warning: replace only with CR 2032 batteries (e.g.: Energizer, Varta, Panasonic).

If there is no blue ray, replace the two AAA batteries immediately, as indicated below:

- remove the lid from the battery compartment;
- remove the old batteries;
- insert the new batteries with the positive (+) and negative (-) sides positioned as indicated by the symbols on the housing;
- replace the battery compartment lid.

Warning: replace only with 1.5 V AAA batteries (e.g.: Duracell, Energizer).

Use of other types of batteries can cause explosions or fires.

Warning: the battery can explode if handled carelessly. Never recharge, remove or position the batteries near sources of heat higher than 100°C (212°F). Do not incinerate.

Dispose of used batteries in compliance with current law.

Maintenance

Pic ThermoKlinik has an infrared sensor located inside the probe. To prevent damage to the sensor, close the thermometer by turning the cover. If there is dust or water inside the probe, clean it delicately with a soft cloth and let it dry before closing.

It is best to occasionally clean the outside of the probe (without touching the sensor) using a 70% alcohol solution.


Switching between °F/°C

While the thermometer is off, press the “M” button (fig. 1 - no. 4) to switch from the “C” scale to “F” and vice versa.

Troubleshooting

Display	Meaning	Action
	There is a general thermometer function problem.	Reset the thermometer by pressing the “ON” and “M” buttons simultaneously and keeping them pressed for 3 seconds. If the thermometer still does not function, contact technical services.
	The operating temperature is outside the thermometer's range of 16°C-36°C (60.8°F - 96.8°F).	Place the thermometer in a room where the temperature is between 16°C and 36°C (60.8°F - 96.8°F). Wait 30 minutes before taking another temperature reading.
	The temperature measured is outside the measurement range 34.4°C - 42.6°C (93.4°F - 108.7°F).	Repeat the measurement and wait for the beep before moving the thermometer.
	The battery is low. Only a few more measurements can be taken.	Replace the old battery with a new one.

Troubleshooting

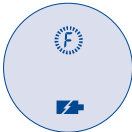Display	Meaning	Action
	The battery is dead. No more measurements can be taken.	Replace the old battery with a new one immediately.
The thermometer does not function or does not respond. The display does not go on or remains on steady.	The thermometer needs to be reset or the battery replaced.	Reset the thermometer by pressing the “ON” and “M” buttons simultaneously and keeping them pressed for 3 seconds or replace the old battery with a new one. If the thermometer still does not function, contact technical services.

Technical data

Measurement interval: 34.4°C - 42.6°C (93.4°F - 108.7°F)

Operating temperature: 16.0°C - 36.0°C (60.8°F - 96.8°F)

Laboratory precision: $\pm 0.2^{\circ}\text{C}$ ($\pm 0.36^{\circ}\text{F}$)

Batteries: 1 3V CR 2032 lithium battery (for the measurement functions)

2 1.5V AAA batteries (only for the light ray)

Classification according to EN60601-1 safety standards:

1. Internally powered equipment
2. Type BF
3. Continuous operation

Warning: refer to the instructions before use and keep them in a safe place.

The thermometer must be kept at room temperature in a dry place, away from any sources of heat and out of direct sunlight. Strong electromagnetic fields could interfere with correct functioning. Keep out of the reach of children. Never leave children alone when taking a temperature reading. If there is a rise in temperature, consult your doctor. Compliant with EEC directive 93/42.

2-year warranty

Artsana guarantees Pic Thermoklinik from material and manufacturing defects for a period of two years from the date of purchase. The warranty applies when the unit is employed for normal home use, in compliance with what is outlined in this user's manual. In no case can Artsana be held responsible for any particular damages, inherent, direct or consequent to the purchase or use of this product or for costs higher than the original cost of the product.

The warranty does not cover any damages due to the batteries, accidents or improper use.

Additional warranties may vary from country to country.

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Artsana S.p.A. – Grandate (CO) Italy